


Name/REF:	TF304, TF304.1, TF304.2, TF304.3, TF304.4
Item No.:	150021, 150028, 150038, 150037, 150040
Intended Use:	Collection of whole blood for preparation of serum eye drops
Classification:	TF304.2: class IIa / rule 2 TF304, TF304.1, TF304.3, TF304.4: class IIa / rule 6 according to annex IX of directive 93/42/EWG
UMDNS:	12-500, medical bags
Description:	Transfer bags TF 304ff provide a blood collection tube and the combined empty bags. The connection tubings are equipped with slide clamps. The products provide either a 17G size needle for punctation or a luer lock site for connection an individually chosen cannula. All configurations provide a predonation-sampling site with adapters of either Holdex or Sarstedt type as well as a site for connecting a container with volume substitution media for the patient/donor.
Drawing:	example TF304.3

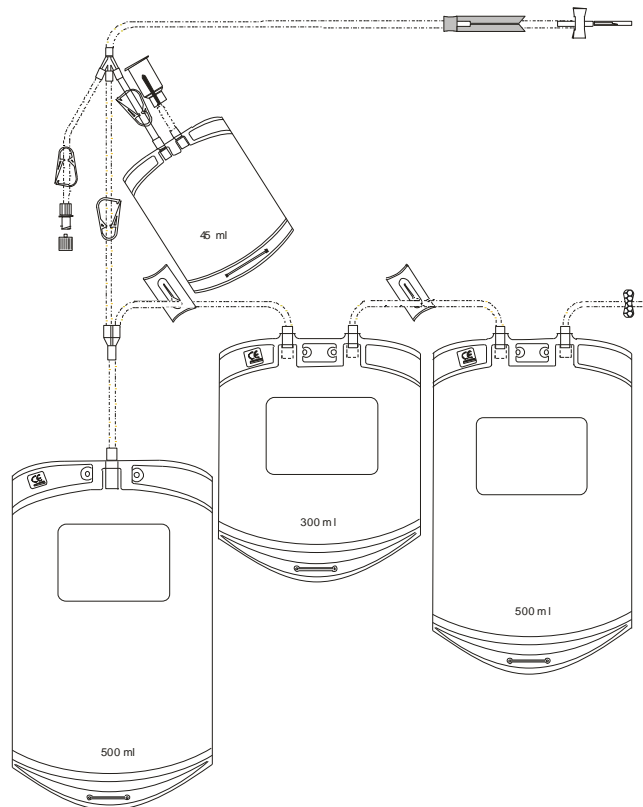

Product configurations:

Features		TF304	TF304.1	TF304.2	TF304.3	TF304.4
Blood withdrawal	17G needle with needle protection	✓	✓	-	✓	✓
	Luer-Lock for connecting a needle	-	-	✓	-	-
Sampling	Sampling bag 45 mL	✓	✓	✓	✓	✓
	Sampling adapter Holdex®	✓	-	✓	✓	-
	Sampling adapter Sarstedt®	-	✓	-	-	✓
Luer-Lock for connecting containers will volume substitution for patient		✓	✓	✓	✓	✓
Collection bag	500 mL	✓	✓	✓	✓	✓
Primary serum bag	300 mL	✓	✓	✓	✓	✓
Final serum bag	300 mL	✓	✓	✓	-	-
	500 mL	-	-	-	✓	✓
	connection site for dilution	-	-	-	✓	✓

Components:

Description	Material
tubing 3,0 x 4,3 mm	PVC, TOTM, Shore 70
Tubing foil	PVC, DEHP, 0,4 mm
needle	17G, 25 mm, backeye, needle protection

Sterilization: gamma radiation according to DIN EN ISO 11137

Packaging: 25 items per box

Storage: Keep dark and dry

Shelf Life: 36 months